



**Compact Massager  
Instructions**

# Compact Massager

# Compact Massager

Please read all instructions carefully to familiarise yourself with your new massager before using.

Save these instructions for further reference.

For any further assistance or information call Customer Services on

**T:** 01227 740066

or visit

**W:** [www.wahl.co.uk](http://www.wahl.co.uk)

**E:** [customer.services@wahl.co.uk](mailto:customer.services@wahl.co.uk)

## CONTENTS

<b>GENERAL USER INFORMATION</b>	Page 4
<b>IMPORTANT SAFEGUARDS</b>	Page 5
<b>DESCRIPTION OF PARTS</b>	Page 7
<b>OPERATING INSTRUCTIONS</b>	Page 8
<b>TROUBLE SHOOTING</b>	Page 8
<b>CLEANING &amp; STORAGE</b>	Page 9
<b>MESSAGE ATTACHMENTS</b>	Page 10
<b>MESSAGE TECHNIQUES</b>	Page 11
<b>GUARANTEE</b>	Page 14

# Compact Massager

## GENERAL USER INFORMATION

### Information about using the operation manual

- ▶ Before using the appliance for the first time, you must read and understand the operation manual completely.
- ▶ Consider the operation manual as part of the product and store it in a safe and accessible place.
- ▶ Include this operation manual if passing the appliance on to a third party.
- ▶ In the case of conflicts with national safety specifications or instructions, the national requirements must be followed.

### Explanation of symbols and notes

**DANGER**

Danger of explosion which may result in serious physical injury or death.

**WARNING**

Warning of possible physical injury or a health risk.

**CAUTION**

Information about danger of material damage.



Note with useful information and tips.

- ▶ Action to be taken by owner.
  1. Carry out these actions in the sequence described.
- List

# Compact Massager

Please retain this leaflet for future reference.

**READ ALL INSTRUCTIONS BEFORE USING.**

## **IMPORTANT SAFEGUARDS**

When using this electrical appliance, especially when children are present, basic safety precautions should always be observed including the following.

### **Requirements**

- ▶ This appliance can be used by children from aged 8 years and above and persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- ▶ Use this appliance only for its intended purpose as described in this manual. Do not use attachments not recommended by the manufacturer.
- ▶ Do not use while bathing or in the shower.
- ▶ Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not

## Compact Massager

place in or drop into water or other liquid.




### **Warning! Injuries due to incorrect handling**

- ▶ In case of pregnancy, diabetes or illness, consult your physician before using this massager.
- ▶ Do not use the massager on:
  - a. Sensitive skin, swollen or inflamed areas, in the presence of poor circulation, where skin eruptions are present, or in the presence of unexplained calf or abdominal pain.
  - b. Varicose veins.
  - c. The front of the neck.
  - d. Genital areas.
  - e. Sleeping or unconscious persons.
- ▶ When using a massage attachment be sure to seat the attachment fully onto the stem and keep it on while massaging.
- ▶ Do not use the massager for more than 20 minutes at one time.
- ▶ Do not use the massager in close proximity to loose clothing or jewellery.
- ▶ Keep long hair away from massager while in use.
- ▶ Do not use the massager without attachments.

# Compact Massager

## DESCRIPTION OF PARTS

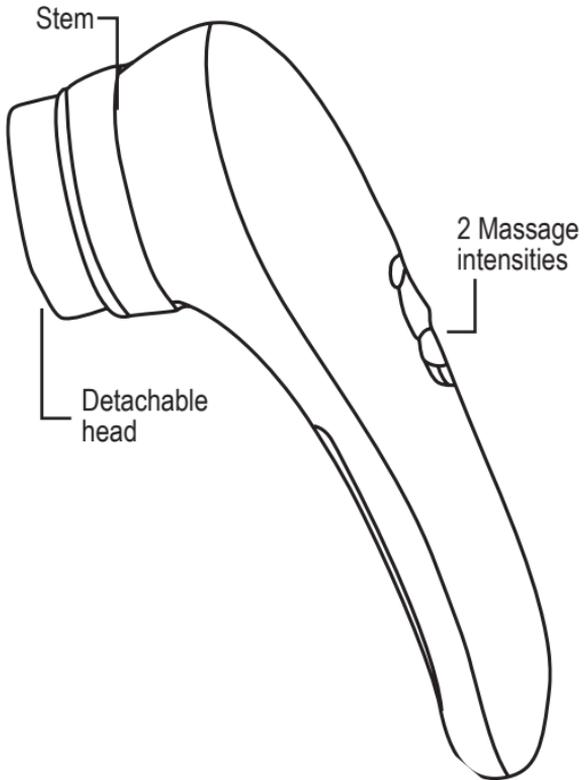

# Compact Massager

## OPERATING INSTRUCTIONS

1. Place the desired massage attachment on the stem of the massager. The attachments easily press on and pull off.
2. Slide the switch to the desired massaging speed -
  - I-Low for soothing massage.
  - II-High for vigorous, deep penetrating action.
3. Apply the massager to the targeted area, using continual movement. A circular motion is often most affective. (See MESSAGE TECHNIQUES.)

## TROUBLESHOOTING

If your massager fails to operate when switched "ON" check your batteries to see if they are properly installed. Positive end should be located as shown in the battery housing. Make sure the switch is in the "ON" position. If it still does not operate, your batteries have probably exceeded their useful life. To maximise battery life, replace both AA batteries at the same time.


## BATTERY REPLACEMENT

When the massager's strength becomes weak, the battery needs to be replaced.

Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (NI-Cad, NiMH, etc) batteries.

1. Press on the "key" and pry off the battery cover (Fig. 1).
2. Remove the batteries and insert two new AA batteries. Make sure the left battery has the positive end up and the right battery has the positive end down (Fig. 2).
3. Replace the battery cover by inserting the bottom of the cover first, then pressing on the upper part.

## Compact Massager

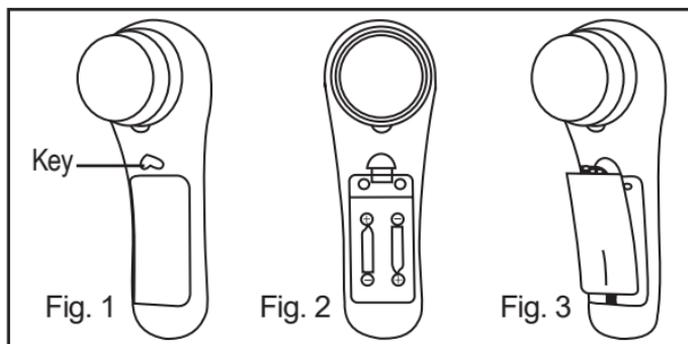

### CLEANING AND STORAGE

**IMPORTANT:** If the unit is going to be stored for a long period of time, be sure to remove the batteries.

After each use, the sponge applicator needs to be cleaned by mild water and allowed to air dry before the next use. If the sponge changes colour or gets mildewed, discard it immediately.

To clean the unit, simply wipe the unit with a dry, soft cloth. Do not immerse in water or use any harsh cleaners or solvents.

### CHANGING ATTACHMENTS

To remove an attachment head, use your thumb and two fingers, and pull away from the massager stem. To insert an attachment head, align the head with the massager stem and press firmly onto the stem.

## Compact Massager

### MASSAGE ATTACHMENTS

\*Some kits may not include all attachments



#### DEEP MUSCLE

Delivers a ripple effect action for large muscle areas.



#### ALL PURPOSE

For general, all body massage.




#### FACIAL SPONGE

Provides beneficial, revitalizing facial massage. (Apply cold cream or lotion to face prior to use. Keep massager away from eyes.)



#### RAISED BUMP

Reaches deep to relieve muscles.

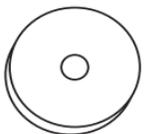

#### ACCUPOINT

Lets you apply deep, pinpoint massage action to the exact body area.

## Compact Massager

### **MASSAGE TECHNIQUES**

#### **SHOULDERS AND NECK —**

Massage each shoulder toward back of neck to point just behind and below ear. Then back to shoulder. Repeat several times.

#### **ARMS —**

Start with a gentle massage of the palms and fingers. Next, using a circular motion, move upward to forearms. Then up to biceps and shoulders.

#### **LEGS —**

If desired, massage bottoms of feet. Then, using a circular motion, massage each leg from calf on up and over thigh.

#### **CHEST —**

Massage from centre of chest outward, and up toward each shoulder.

#### **BACK —**

Massage from top of spine to bottom, and from centre outward.

#### **ABDOMEN —**

Relax abdominal muscles by massaging from centre outward, using a gentle, circular motion.

#### **JOINTS —**

Relieve tired, aching knees and elbows with smooth, circular strokes.

## Compact Massager

### **FACE —**

Gently massage from the nose outward, taking care to avoid contact with the eye.

### **SCALP —**

Start at the back of the neck and massage up to the top of the head.

Using this up and down motion, massage across the scalp.

## Compact Massager

### **BATTERY DISPOSAL**

#### **Battery Removal and Disposal**

Do not dispose of this product or the battery in the normal household waste.

At the end of your products life the battery must be removed and recycled/disposed of correctly and not disposed of with the product.

Batteries should be disposed of at an official collection point in an environmentally safe way or taken to a local collection facility where they will be dealt with safely and sent for the appropriate recycling.

### **WASTE ELECTRICAL & ELECTRONIC EQUIPMENT DIRECTIVE 2002/96/EC (WEEE)**



At the end of this products life it should not be disposed of with general household waste.

Care should be taken to dispose of it in a suitable manner in accordance with your Local

■ Authority regulation.

## Compact Massager

### **GUARANTEE**

This product is guaranteed against defects in material and workmanship for a period of 2 years from the date of original purchase or receipt as a gift.

Should this product become defective during the guarantee period, return it to the store of purchase together with your proof of purchase for repair or replacement

Alternatively within the guarantee period you can return the product to Wahl (UK) Ltd who will repair any such defect or elect to replace the product or any part of it, without charge, provided that there is proof of purchase. Should a replacement be offered this will not extend the original guarantee period.

The guarantee becomes invalid in the case of misuse, alteration or repair by unauthorised persons.

This product is designed for domestic use only and use within a professional environment will invalidate the guarantee.

This guarantee in no way affects your rights under statutory law in the United Kingdom.

## Compact Messenger

WAHL (UK) LTD  
STERLING HOUSE  
CLIPPER CLOSE  
RAMSGATE  
KENT  
CT12 5GG

CUSTOMER SERVICES T: 01227 740066  
[www.wahl.co.uk](http://www.wahl.co.uk)

Form No. ZB343 Rev.2 08/18